

BAC-PACK™ Bougie-Aided Cricothyroidotomy Pack

REV091323



A comprehensive solution that allows the user to perform a surgical cric in any environment

The North American Rescue® Bougie-Aided Cricothyroidotomy Pack - BAC-Pack™ is a comprehensive solution that allows the user to perform a surgical cric in any environment. With the components arranged by order of use and organized with elastic loops, the BAC-Pack is an easy-to-use kit in a small pouch made of durable 500D nylon.

The BAC-Pack™ utilizes a purpose-built tube that has a pre-inserted bougie for rapid application, and the bougie is fully adjustable based on the patient's or provider's needs. The excess bougie can also be conveniently stored in the one-of-a-kind bougie cap/stopper. The tube meets the TCCC recommendations for the open surgical technique (flanged and cuffed airway cannula of less than 10 mm outer diameter, 6-7 mm internal diameter, and 5- 8 cm of intratracheal length).

The BAC-Pack™ also comes with the latest innovation in safety scalpels, the #10 Protected Scalpel, for effortless penetration through the cricothyroid membrane and the NAR® Tracheal Hook to maintain the opening without worry. A single ChloroPrep™ swabstick, an alcohol pad, and a gauze pad are included for skin prep, and a 10cc syringe and tube securing strap are included for finishing the procedure.

Features:

- All-in-one cricothyroidotomy pack for use in any environment
- Pre-inserted bougie in tube allows for a more rapid, easy procedure
- Tube meets TCCC recommendations for bougie-aided open surgical technique
- #10 Protected Scalpel included for safety and ease of use
- Small, durable pouch keeps the items organized and ready to use

Specifications:

- Packaged:
H 8 in. x W 4 in. x D 1.75 in.
- Weight: 5.2 oz

BAC-Pack™

ITEM #	NSN #
10-0052	----

SINGLE USE ONLY
 Rx ONLY
 NOT MADE WITH NATURAL RUBBER LATEX

Nylon Color COY



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277

PRODUCTS
 WITH A **MISSION®**


BAC-PACK™ Bougie-Aided Cricothyroidotomy Pack

REV091323



Kit Contents:

- 1x Rolled Nylon Pouch
- 1x Alcohol Prep Pad
- 1x Chloraprep Swabstick
- 1x 2-Pack Gauze Pad (4 in. x 4 in.)
- 1x #10 Protected Scalpel
- 1x NAR® Tracheal Hook
- 1x NAR® Tracheostomy Kit with Bougie Introducer (Cuffed, 6mm)
- 1x Tracheal Tube Securing Strap
- 1x 10cc IV Syringe

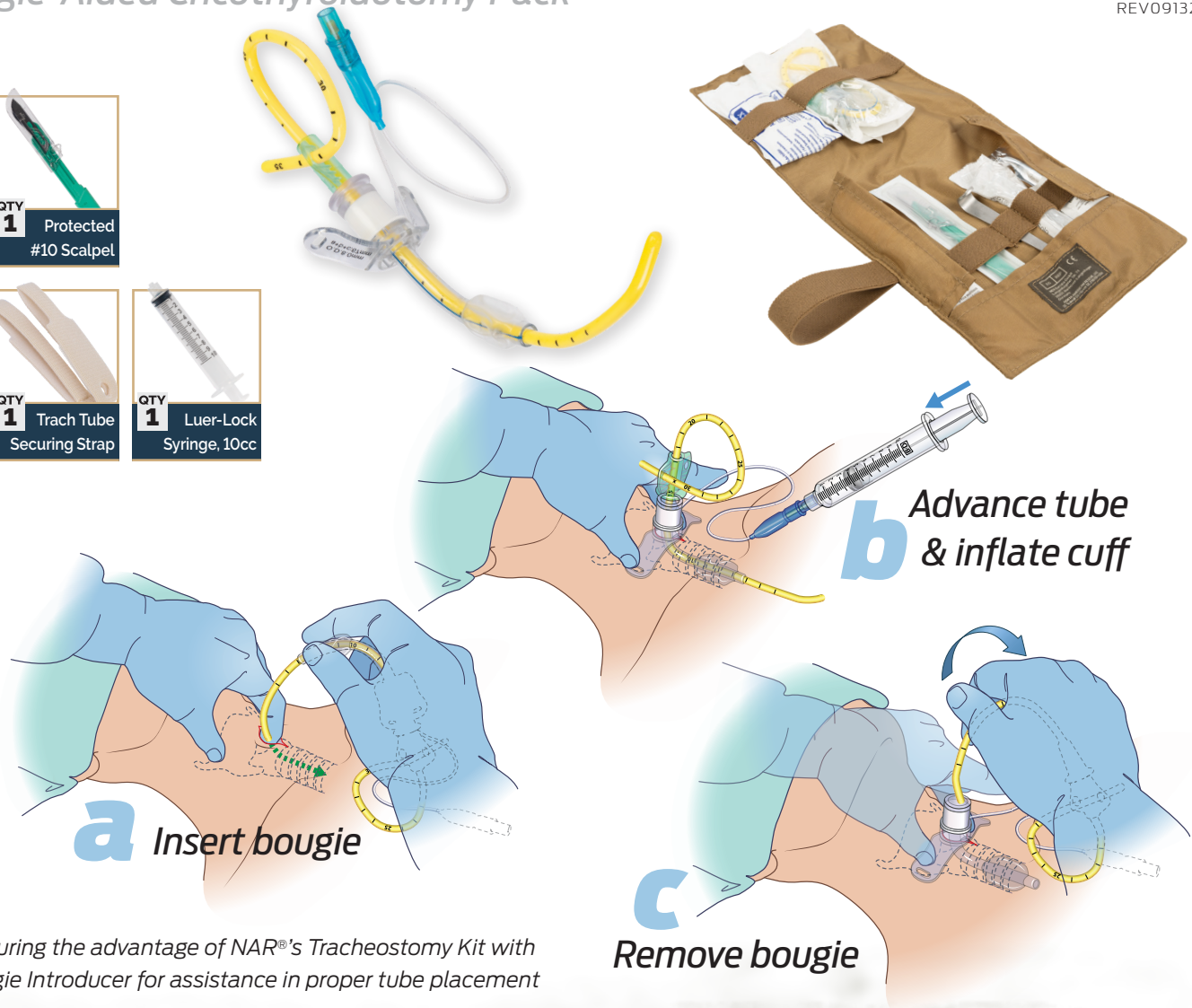

featuring the advantage of NAR's Tracheostomy Kit with Bougie Introducer for assistance in proper tube placement



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277

PRODUCTS
WITH A **MISSION®**