

SAM
MEDICAL

SAM[®] XT

Extremity Tourniquet

**CLICK.
TWIST.
SECURE.**

PROVEN TO SURVIVE

Engineered for rapid application, SAM XT's innovative design enables easier, faster training and intuitive use.



FEATURED TECHNOLOGY

TRUFORCE™ BUCKLE

Slack is the main cause of failed tourniquet application. The TRUFORCE™ Buckle technology auto-locks to eliminate nearly all tourniquet slack by incorporating innovative baseline force control, activating the locking prongs.






SAM[®]XT


ENGINEERED TO PRESERVE LIFE

- 1 TIME BAND**
Weather-resistant hook and loop system allows for easy recording time of application.
- 2 WINDLASS C-HOOKS**
Secures the Windlass Rod using an intuitive angular design, making windlass locking easier and more efficient.
- 3 TRUFLEX™ PLATE**
Durable nylon construction distributes force as targeted pressure for maximum strength and flexibility, able to bend on smaller extremities without breaking.
- 4 TRUFORCE™ BUCKLE**
Innovative auto-lock technology activates locking prongs, eliminating slack typically found in other tourniquet applications.

- 5 HIGH IMPACT POLYMER WINDLASS**
Injection molded from strong, lightweight, heat and UV stable High Impact Polymer.
- 6 PRECISION STRAP**
Velcro-dependent tourniquets can fail when wet or dirty. The all-weather Precision Strap is made of single-layer nylon webbing. It features intervalled holes that lock with the TRUFORCE™ Buckle to create a consistent, efficient primary security system. The secondary hook and loop system secures the remaining strap, adding stability and safety.
- 7 PRECISION TAB**
Intuitive design allows for quicker, easier tourniquet application. Material is highly visible in low-light situations.

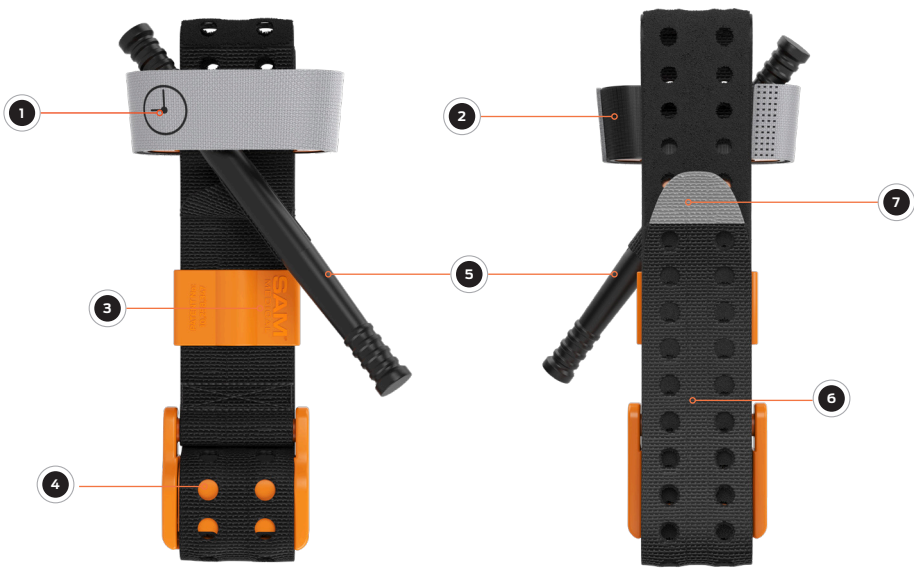

MODELS

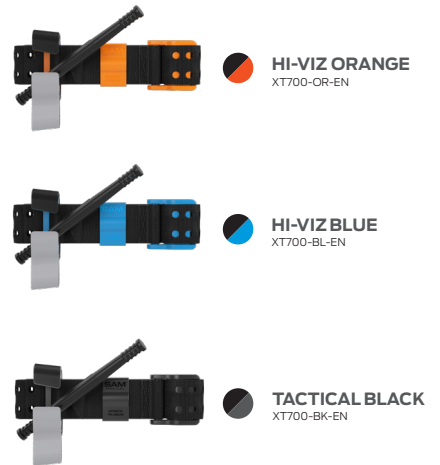

PRODUCT SPECIFICATIONS

PACKAGED SIZE & WEIGHT*

Length: 2.5 in (6.3 cm)
 Width: 1.5 in (3.8 cm)
 Height: 7.25 in (18.4 cm)
 Weight: 3.4 oz (96.4 g)

OPERATING SPECIFICATIONS

TRUFORCE™ Buckle technology auto-locks to engage strap holes within circumferential range:
 5.75 in - 32.5 in (14.6 cm - 82.6 cm)

MODELS AVAILABLE

Tactical Black (XT700-BK-EN)
 HI-VIZ Orange (XT700-OR-EN)
 HI-VIZ Blue (XT700-BL-EN)

GLOBAL SALES

+1 503.639.5474
sammedical.com

SAM Medical Products®
 12200 SW Tualatin Rd., Ste 200
 Tualatin, OR 97062 USA

U.S. Patent No. 10,258,347 | 10,363,046
 *Size and weight vary by configuration and manufacturing process.
 SAM® is a registered trademark in the United States and other countries.
 ©2024 SAM Medical Products. All rights reserved.

DISTRIBUTED BY

Life-Assist
 Helping Heroes Save Lives