



# STAIR CHAIR 4BELL SILVER+

## FEATURES

- **Frame in aluminum** makes the chair lightweight and resistant at the same time
- **Two telescopic handles at the front and two foldable ones behind the backrest** help lifting and transport operations
- **Rear wheels' locking mechanism** makes sure the chair does not slide on a flat surface
- **Semirigid seat and soft backrest** are built to increase patient comfort
- **Foldable footrest** increases patient stability
- **Easy to clean PVC seat and backrest** are resistant to chemicals
- **Pivoting front wheels** makes the chair more maneuverable.



Available models	<b>4BELL SILVER+</b>
Frame color	
Seat color	
Code	<b>ST10500</b>

**4BELL SILVER+** is a stair chair designed to transfer or evacuate patients with mobility limitations from multi-storey buildings. The chair is FDA compliant and built in accordance to the specifications set by EN 1865-1.

If the stair chair is installed on the 4Bell Max+ bracket, it will comply with 10G crash test according EN 1789 standard.

The frame of the chair is made in aluminum for increased resistance. Its total weight is 9,4 kg (20.7 lb).

397 lb / 180 kg



Load Capacity

20.7 lb / 9,4 kg



Lightweight



Foldable Chair




Easy to Store

## TECHNICAL DATA

Length with handles retracted	76 cm / 29.9 in
Length with handles extracted	112 cm / 44 in
Width	49 cm / 19.8 in
Height	90,5 cm / 35.6 in
Dimensions when folded	49 x 90,5 x 29 cm / 19,8 x 35,6 x 10.5 in
Weight	9,4 kg / 20.7 lb
Load capacity	180 kg / 397 lb
Front wheels	Ø 10 cm / 3.93 in
Rear wheels	Ø 15 cm / 5.90 in
Materials	Al, Steel, PU, PVC

### THE PRODUCT INCLUDES

Restraint belts.

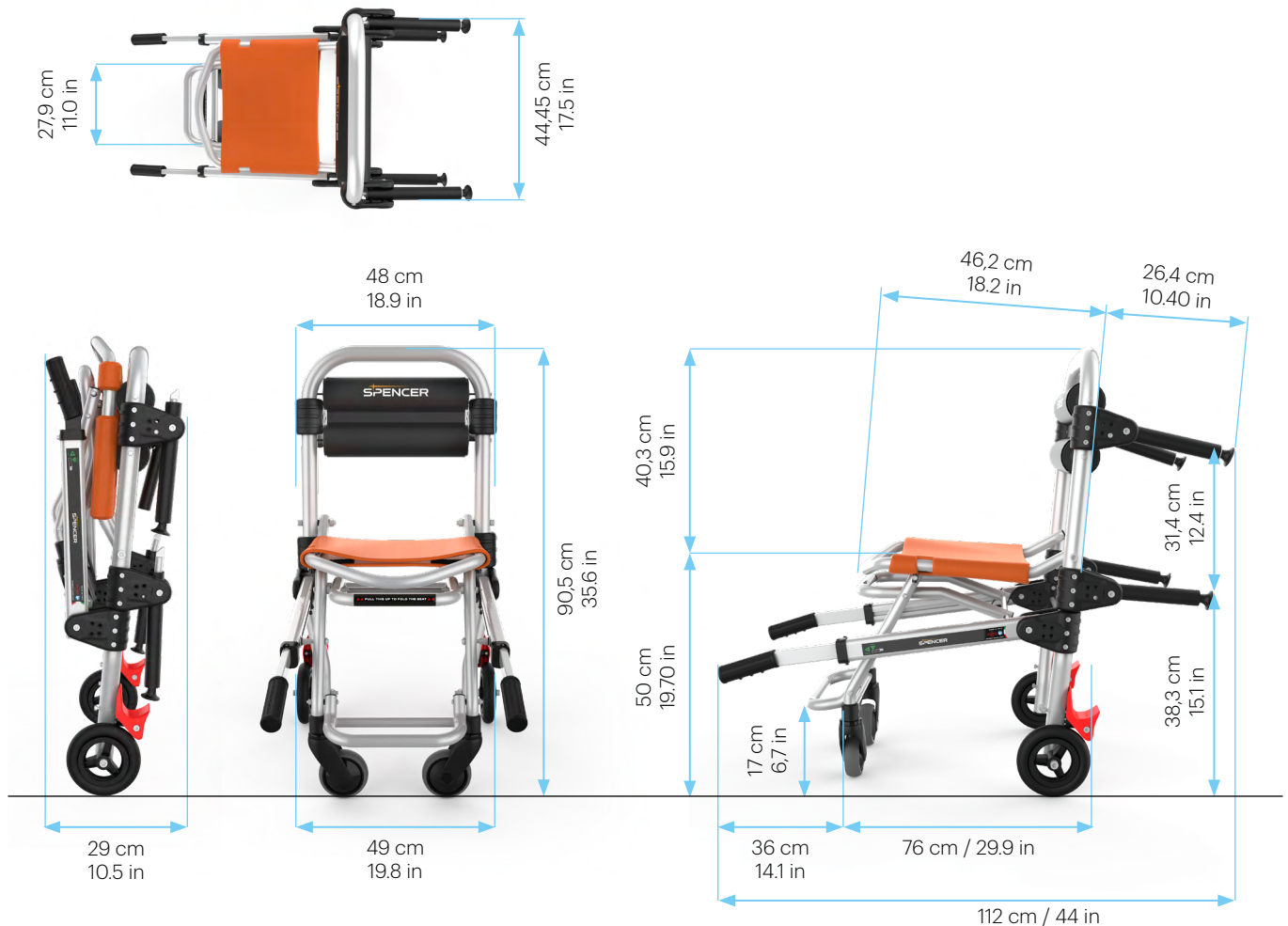



## ACCESSORIES



ST10453

4BELL+ MAX  
10G tested Fastening system

- It is installed in ambulances, special vehicles or airplanes
- It helps securing the 4Bell Silver+ in foldable position to either the wall or the floor
- When 4Bell Silver+ is installed on the 4Bell Max+ bracket, it will comply with 10G crash test according EN 1789 standard.



ST10452	XL Seat
ST10450	Armrests
ST10451	Rear auxiliary handles

Rev.0 (14/03/2022)

©Copyright - This document is property of SPENCER ITALIA S.r.l and can't be reproduced, used and disclosed without the written approval of SPENCER ITALIA SRL.