

SAM
MEDICAL

SLISHMAN Traction Splint (STS Gen 2 / STS-C)

**QUICK.
EASY.
VERSATILE.**



TRACTION FROM THE HIP

STS GENERATION 2 / STS-C

- One size fits all, from small children to large adults
- Lightweight and compact
- Can be applied in under a minute with minimal training
- No assembly required
- No extension beyond the foot
- Maximizes patient comfort with well padded pressure points
- Radiolucent
- Connected components prevent lost parts
- Can still be applied on patients with lower leg trauma



TRACTION FROM THE HIP

1 NEOPRENE OUTER TUBE STRAP

Hook & loop strap to be secured at the ankle or proximal calf just distal to the patella. Both STS Gen 2 and STS-C fits pediatric and adult patients. STS-C is proximal calf placement only.

2 PADDED GROIN STRAP

Secures at the patient's hip. Padded for patient comfort.

3 TRACTION CORD


Red-handled traction cord applies traction to patient. Traction is locked by moving cord into V-lock on red traction pole.

4 PRESSURE WRAP

Multi-function wrap can be used to buddy splint patient's legs or feet, apply bleeding control pressure to wounds or as a wrap for securing other splints to patient.

5 CARRY CASE

STS label for easy identification of splint. Integrated daisy chain for securing to packs or bags.

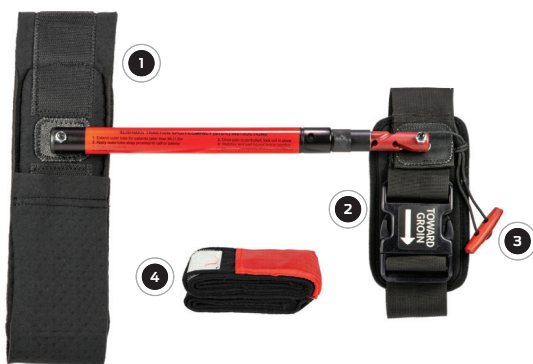

STS GEN 2 (SLISHMAN TRACTION SPLINT GENERATION 2)

Advanced design and versatile features make the STS Gen 2 a must-have device for any responders treating femur fractures. Extended length allows for no extension beyond the foot to interfere in patient transport.



CARRY CASE WITH PALS LOOPS



TRAINING VIDEO

STS-C (SLISHMAN TRACTION SPLINT COMPACT)

Extremely compact and lightweight, the STS-C is designed for situations where medical equipment must be packed into difficult or remote locations. Ideally suited for search and rescue, special operations, military medics, wildland firefighters, backcountry skiers, or anyone trying to minimize pack weight and cube space.



CARRY CASE WITH MOLLE SNAPS



TRAINING VIDEO

PRODUCT SPECIFICATIONS

PACKAGED SIZE & WEIGHT (Dimensions in Bag)

STS Gen 2

Height: 22 in (55.9 cm)
Width: 3 in (7.6 cm)
Depth: 3 in (7.6 cm)
Weight: 18 oz (0.51 kg)

Patents: 10,517,750 & 11,324,624
CE marked / NSN pending

STS-C

Height: 14 in (35.6 cm)
Width: 3 in (7.6 cm)
Depth: 3 in (7.6 cm)
Weight: 15.5 oz (0.44 kg)

CE marked
US, EP & Canadian patents
NSN: #6515-01-694-7603

MODELS AVAILABLE

STS Gen 2 (10-5000)
STS-C (10-1270)

GLOBAL SALES

503.639.5474

Distributed By

SAM Medical
12200 SW Tualatin Rd, Ste 200
Tualatin, OR 97062 USA
sammedical.com